

Validation of Snort-spit Saliva in Detecting COVID-19 Using RT-PCR and Rapid Antigen Detection Test

Ryner Jose D. Carrillo, MD, MSc (Clin Epid),^{*1,2} Abigail D. Sarmiento, MD,^{*1} Mark Anthony C. Ang, MD, MoS,^{1,2} Michele H. Diwa, MD,^{1,2} Cecile C. Dungog, MD, PhD,¹ Daniel I. Tan, MD,¹ Jan Alexis C. Lacuata, MD,¹ Jacob Ephraim D. Salud, MD,¹ Ramon Antonio B. Lopa, MD,^{1,2} John Mark S. Velasco, MD, MPH, MSc, DTM&H,³ Regina P. Berba, MD, MSc^{1,2} and Charlotte M. Chiong, MD, PhD,^{1,2,3}

¹Philippine General Hospital, University of the Philippines Manila

²College of Medicine, University of the Philippines Manila

³National Institutes of Health, University of the Philippines Manila

ABSTRACT

Objective. To determine the diagnostic accuracy of self-collected snorted and spit saliva in detecting COVID-19 using RT-PCR (ssRT-PCR) and lateral flow antigen test (ssLFA) versus nasopharyngeal swab RT-PCR (npRT-PCR).

Methods. One hundred ninety-seven symptomatic subjects for COVID-19 testing in a tertiary hospital underwent snort-spit saliva self-collection for RT-PCR and antigen testing and nasopharyngeal swab for RT-PCR as reference. Positivity rates, agreement, sensitivity, specificity, and likelihood ratios were estimated.

Results. Estimated prevalence of COVID-19 using npRT-PCR was 9% (exact 95% CI of 5.5% - 14.1%). A higher positivity rate of 13% in the ssRT-PCR assay suggested possible higher viral RNA in the snort-spit samples. There was 92.9% agreement between ssRT-PCR and npRT-PCR (exact 95% CI of 88.4% to 96.1%; Cohen's Kappa of 0.6435). If npRT-PCR will be assumed as reference standard, the estimated Sensitivity was 83.3% (exact 95% CI of 60.8% to 94.2%), Specificity 93.9% (exact 95% CI of 89.3% to 96.5%), Positive predictive value of 57.7% (exact 95% CI of 38.9% to 74.5%), Negative predictive value of 98.2% (exact 95% CI of 95% to 99.4%), positive likelihood ratio of 3.65 (95% CI of 7.37 to 24.9), negative likelihood ratio of 0.178 (95% CI of 0.063 to 0.499). There was 84.84% agreement (95% exact CI of 79.1% to 89.5%; Cohen's Kappa of 0.2356) between ssLFAvs npRT-PCR, sensitivity of 38.9% (exact 95% CI of 20.3% to 61.4%), specificity of 89.4% (exact 95% CI of 84.1% to 93.1%), PPV of 26.9% (95% CI of 13.7% to 46.1%), NPV of 93.6% (exact 95% CI of 88.8% to 96.4%), LR+ of 3.67 (95% CI of 1.79 - 7.51), LR - of 0.68 (95% CI of 0.47 - 0.99).

Conclusion. Our data showed that snort-spit saliva RT-PCR testing had acceptable diagnostic performance characteristics and can potentially be used as an alternative to the standard nasopharyngeal/oropharyngeal swab RT-PCR test for COVID-19 in certain situations. However, our data also showed that snort-spit saliva antigen testing using lateral flow assay did not offer acceptable performance.

Key Words: SARS-CoV-2, reverse transcription polymerase chain reaction, saliva

INTRODUCTION

In the era of COVID-19 or SARS-CoV-2 pandemic, determining infection and infectivity by viral RNA reverse transcription polymerase chain reaction (RT-PCR) has been the diagnostic method of choice. Studies have tested different samples- nasopharyngeal (NP) or oropharyngeal (OP) swabs, bronchial lavage, whole saliva, and deep throat saliva. The nasopharyngeal swab RT-PCR has been considered by the World Health Organization as the sample type of choice for diagnostics¹. However, the collection of NP specimens can be regarded as invasive and causes discomfort

**Dr. Carrillo and Dr. Sarmiento share primary authorship for this manuscript.*

Corresponding author: Ryner Jose D. Carrillo, MD, MSc (Clin Epid)
University of the Philippines-Philippine General Hospital
College of Medicine
University of the Philippines-Manila
Taft Avenue, Ermita Manila 1000, Philippines
Email: rdcarrillo@up.edu.ph

resulting in coughing in most patients and increasing the risk of nosocomial spread of respiratory viruses. The process of obtaining these samples exposes the healthcare professional to infection such that the process will require a high level of personal protective equipment (PPE).

Salivary glands have been proven to be reservoirs for COVID-19 asymptomatic infection.² Saliva or oral fluid can easily be obtained by asking the patient to cough out saliva from the throat to container sample for RT-PCR.³ Serial self-collected saliva testing yielded positive in 91.7% (11/12) patients showing that a self-collected saliva is a viable option in detecting 2019-nCoV.⁴ In a 2020 study by Wyllie, saliva and nasopharyngeal samples were compared, and they detected SARS-CoV-2 in saliva from two subjects who were negative by nasopharyngeal swabs.⁵ They also found that saliva yielded higher detection sensitivity and consistency throughout infection. A meta-analysis by Czumbel in 2020 found that sensitivity of saliva tests is high at 91% (CI 80–99%) while sensitivity for nasopharyngeal swab is at 98% (CI 89–100%) among previously confirmed COVID-19 patients.⁶

An affordable, easily deployable, rapid, and accurate self-collected test for COVID-19 may be a game-changer in the fight against this pandemic. Rapid antigen tests using lateral flow immunoassays have been explored as a possible alternative in settings where RT-PCR may not be readily available. One operational drawback is that these tests rely mostly on nasopharyngeal specimens. Thus, professional sampling and the use of personal protective equipment during sampling and processing are still necessary. Using more readily available and easily collected samples such as saliva may help decrease the need for protective equipment and decrease the risk of exposure to healthcare workers. Saliva being used in antigen tests may also be used in large-scale and point-of-care settings.⁷

The BIOCREREDIT COVID-19 Ag test (RapiGEN Inc., Gyeonggi-do, 14119, Korea) is a rapid antigen test already being used in the community. It has a published product description of having 94.4% sensitivity and 100% specificity in a population of 134 and a disease prevalence of 26.5%.⁸ Table 1 shows the pooled results from European and South American subjects.

This study aimed to validate self-collected snorted and spit saliva in detecting COVID-19 using RT-PCR (ssRT-PCR) and lateral flow antigen test (ssLFA) versus nasopharyngeal swab RT-PCR (npRT-PCR).

Table 1. Biocredit Rapid antigen test product insert results

	+ RT PCR	- RT PCR	Total
Biocredit COVID Ag +	34	0	34
Biocredit COVID Ag-	2	134	136
Total	36	134	136

METHODS

This was a cross-sectional diagnostic accuracy study, with prospective recruitment, using the STARD classification⁹. This research was conducted from September to November 2020 at the Philippine General Hospital (PGH). The study secured approval from the Institutional Biosafety and Biosecurity Committee and the University of the Philippines Manila Research Ethics Board before commencing (Registration: UPMREB 2020-301-01). Data collection was planned before the index test, and the reference standard was performed.

Two hundred eight (208) symptomatic adults (19-65 years old) suspected and probable cases of SARS-CoV-2 infections, using the WHO case definitions, were enrolled in the study.¹ These included health care workers and individuals who consulted at the outpatient swabbing station of the hospital. All were ambulatory. Those who were severely symptomatic, non-ambulatory, and requiring O₂ support or intubation were excluded from the study.

Specimen collection and transport

A designated health care worker (HCW) supervised the snort-spit saliva sample collection while trained medical technologists obtained the nasopharyngeal swab sample according to standard operating procedures. The HCW instructed subjects to snort and then spit 1-2 mL of the whole saliva directly into two different sterile tubes: the first tube contains viral transport medium (VTM) to be processed for RT-PCR testing, and the second tube contains the assay diluent for the lateral flow immunoassay.

Following the snort-spit saliva collection, the medical technologist performed the nasopharyngeal swabbing on bilateral nostrils using one pledget. The pledget was then placed in the tube with VTM. Triple-packed VTMs of both snort-spit saliva and NP swabs were then immediately transported to the PGH Molecular Laboratory for processing.

Rapid antigen testing of saliva specimens

BIOCREREDIT COVID-19 Ag test (RapiGEN Inc., Gyeonggi-do, 14119, Korea) was used for the rapid testing part. It was a lateral flow immunochromatographic assay intended for the qualitative detection of SARS-CoV-2 antigen.⁸ This assay was designed for use with nasopharyngeal swab samples but was used in this study on saliva samples. The HCW performed the rapid test immediately after saliva collection. Three to four drops of the obtained snort-spit saliva and assay diluent mixture were squeezed into the test device's sample well. Results were considered positive if two bands appear, one red control line and one black test line, after a standard waiting time and regardless of color intensity.

Specimen processing for RT-PCR

The NP swab and snort-spit saliva samples were aliquoted and heat-inactivated at 65 Celsius for 10 minutes or chemically inactivated using the lysis buffer before RNA extraction. SARS-CoV-2 viral RNA amplification and detection were performed following standard operating procedures. Viral RNA extraction was done using GenAmplify™ Viral RNA Purification Kit (The Manila HealthTek, Inc.). SARS-CoV-2 RT-PCR detection was then performed using either Novel Coronavirus (2019-nCoV) Nucleic Acid Diagnostic Kit (Sansure Biotech Inc.) or Real-Time Fluorescent RT-PCR Kit for Detecting SARS-CoV-2 (BGI Genomics Co. Ltd.). Given an NP swab and snort-saliva specimen pair, it was ensured that the same SARS-CoV-2 RT-PCR detection kit was used. All kits used were validated before use.

RESULTS

A total of 208 subjects were recruited and were able to give their informed consent to participate. However, two participants got swabbed first before snort-spit saliva collection, one participant had her NP RT-PCR processed using the automated process, and eight had invalid RT-PCR results. A total of 11 participants were excluded from the final set, and the final number was 197. Participants' median age was 32 years (range 22–64), and 123 participants (62.4%) were female. The prevalence of COVID-19 infection in the study cohort, based on the reference method of npRT-PCR, was 9% (18/197, 95% CI of 5.5% – 14.1%). The snort-spit saliva RT PCR test (ssRT-PCR) and snort-spit antigen test (ssLFA), however, showed higher positivity rates, both 13.1% (26/197), compared to npRT-PCR.

The ssRT-PCR showed a substantial agreement of 92.89% (95% exact CI of 88.4% to 96.1%; Cohen's Kappa of 0.6435) with npRT-PCR (see Table 2). The ssRT-PCR has a positive likelihood ratio of 13.66 and a negative likelihood ratio of 0.178. This means that if the pretest probability of a case is 9%, and the ssRT-PCR result is positive, the post-

test probability will be 57.7%. If the pre-test probability is 9%, and the ssRT-PCR result is negative, the post-test probability will be 1.7%. The sensitivity is 83.3% (95% exact CI of 60.8% to 94.2%) and the specificity is 93.9% (95% exact CI of 89.3% to 96.5%).

The ssLFA showed an only fair agreement of 84.77% (95% exact CI of 79.1% to 89.5%; Cohen's Kappa of 0.2356) with npRT-PCR (see Table 3). The ssLFA has a positive likelihood ratio of 3.67 and a negative likelihood ratio of 0.68. This means that if the pretest probability of a case is 9%, and the ssLFA result is positive, the post-test probability will be 26.6%. If the pre-test probability is 9%, and the ssLFA result is negative, the post-test probability will be 6.3%.

Further analysis of the data (Figure 1) shows considerable clustering of ssRT-PCR cycle threshold (Ct-values) towards the higher Ct-values in samples read as negative by the reference standard (npRT-PCR), suggesting some form of systematic error.

Comparison of the distribution of npRT-PCR and ssRT-PCR Ct-values with positive and negative ssLFA results shows that ssLFA misclassified a considerable proportion of cases with detectable viral RNA by RT-PCR in both the npRT-PCR and ssRT-PCR and that the misclassification error is random as it does not demonstrate any pattern across the range of Ct-values (number of PCR cycles to achieve the desired *signal*). (Figure 2 A and B.)

DISCUSSION

Nasopharyngeal swab RT PCR is currently considered the diagnostic method of choice for SARS COV 2 testing by the WHO guidelines.¹ Due to its invasive technique, high risk for coughing, increased risk of exposure to HCWs, and higher demand for protective equipment for sample

Table 2. Snort-spit saliva RT-PCR vs. Nasopharyngeal swab RT-PCR for SARS-Cov-2

	npRT-PCR +	npRT-PCR -	Total
ssRT-PCR +	15	11	26
ssRT-PCR -	3	168	171
Total	18	179	197

Table3. Snort-spit saliva antigen test vs. Nasopharyngeal swab RT-PCR for SARS-Cov-2

	npRT-PCR +	npRT-PCR -	Total
ssLFA +	7	19	26
ssLFA -	11	160	171
Total	18	170	197

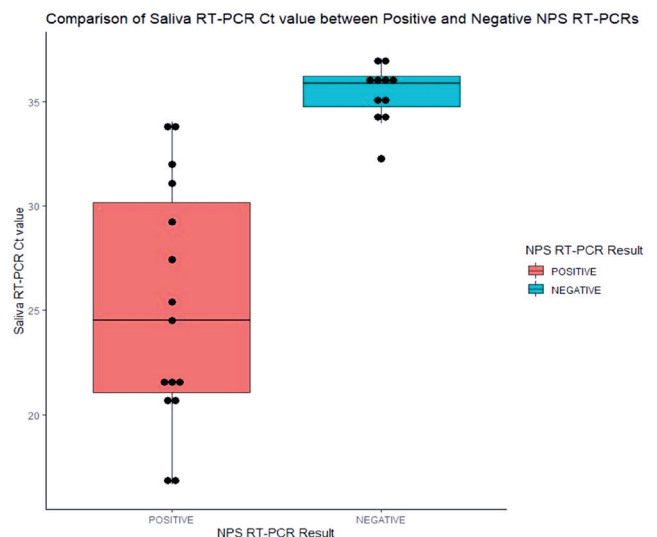


Figure 1. Comparison of Saliva RT-PCR Ct values between Positive and Negative NPS RT-PCRs.

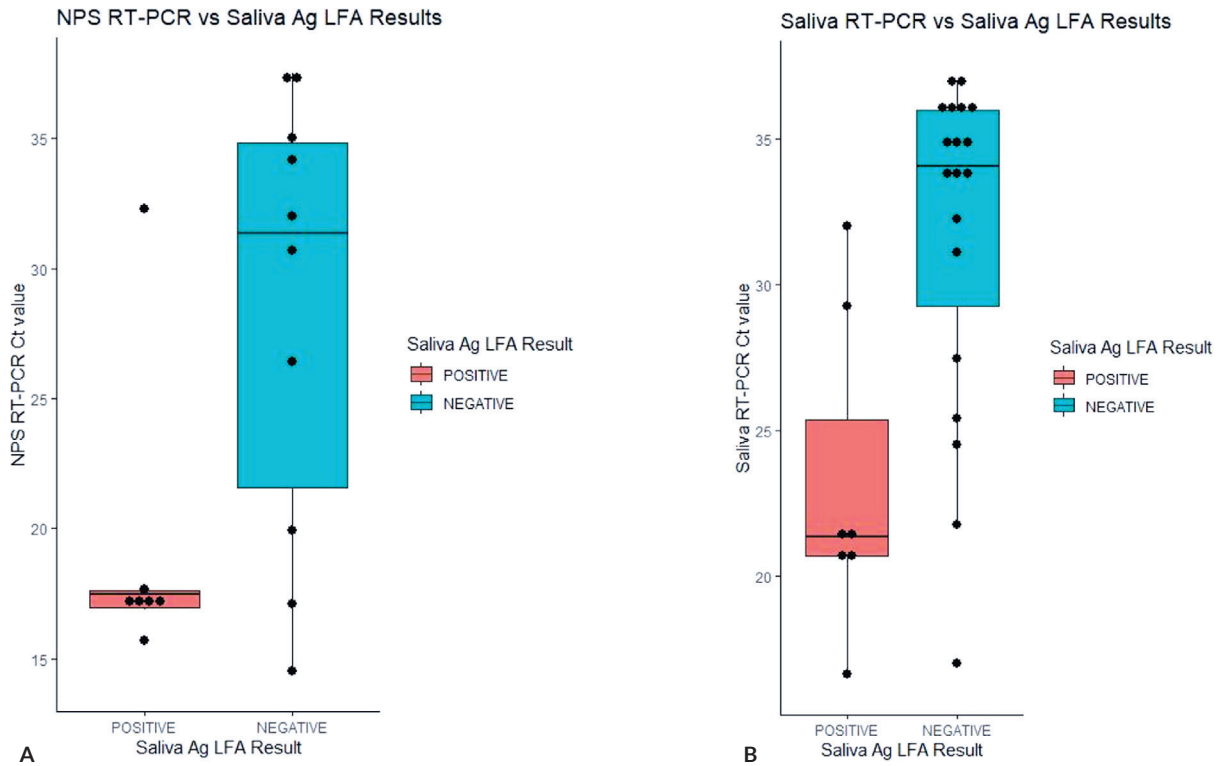


Figure 2. Comparison of snort-spit saliva rapid antigen testing (ssLFA) with (A) snort-spit saliva RT-PCR (ssRT-PCR) and (B) nasopharyngeal RT-PCR (npRT-PCR).

collection, several studies have explored saliva- or sputum-based testing. The current study's difference is that it utilizes the snort-spit technique of collection, aiming to sample both the nasopharyngeal mucous and whole saliva.

Snort-spit Whole Saliva RT-PCR

Using the npRT PCR as the reference standard, the diagnostic sensitivity and specificity of ssRT-PCR in detecting COVID19 infection is high at 83.3% and 93.9%, respectively. This is consistent with previous cross-sectional studies by Pasomsab (2020), with saliva having sensitivity and specificity of 84.2% and 98.9%, respectively, and a more recent study by Altawalsh (2020) showing sensitivity and specificity of 83.43 % (95 % CI: 79.07–87.20) and 96.71 % (95 % CI: 94.85–98.04 %), respectively.^{10,11} However, the presence of a wide variation in the assay performance based on several more studies results, the WHO does not yet recommend using saliva as the sole specimen for routine clinical diagnostics.¹

The positive and negative likelihood ratio estimates in this study for ssRT-PCR suggest usefulness in an actual clinical setting. For example, in a patient assigned a pretest probability of 50%, that is, the clinical decision-maker is at an equipoise, a positive ssRT-PCR can increase the post-test probability to 93.2%, and a negative ssRT-PCR can decrease the post-test probability to 15.1%, enough to be used as the basis for clinical decisions.

The current study provides some evidence that the detection rate of SARS-COV-2 could be higher in snort-spit samples than nasopharyngeal swabs (13.3% vs. 9% positivity rate, respectively). If cross-contamination at any point in the processing can be ruled out, a higher positive test rate will suggest an actual increase in detection rate. Furthermore, because there was considerable clustering of the Ct-values towards the higher end of the assay range, in the group of positive samples by ssRT-PCR but negative by npRT-PCR, the npRT-PCR assay could be systematically misclassifying those sample with very low viral copies. Thus it is very tempting to consider the possibility that ssRT-PCR could, in reality, be more sensitive than npRT-PCR in SARS-CoV-2 detection. If indeed this could be true, it would be easy to see that such will be the case as samples collected by snorting and spitting collects fluid from a larger area, including the nasal, nasopharyngeal, oropharyngeal, and oral cavities, and the total fluid collected is far larger in volume compared to nasopharyngeal and oropharyngeal swabbing. Further investigation using carefully designed experiments is suggested to explore this assertion. The pre-amplification extraction method and validation of SARS-CoV-2 RT PCR interpretation criteria specifically for snort-spit samples may be pursued as a future direction to provide additional useful information to strengthen the evidence for the potential of snort-spit as an alternative sample for SARS-CoV-2 detection.

Snort-spit Saliva Rapid Antigen Assay (Lateral Flow Assay)

We explored the possibility of snorting and spitting mucus and saliva to a collection bottle and testing it on a lateral flow device for rapid antigen detection using the BIOCREDIT COVID-19 Ag test (RapiGEN Inc., Gyeonggi-do, 14119, Korea). With a snort-spit sample, there was only fair agreement and poor sensitivity and specificity when using antigen than what is published in the product insert, which uses nasopharyngeal swab. The current study employed a modified sample processing method to perform the antigen test, which may be a limitation. Hence, the utility of snort-spit saliva samples using rapid antigen detection may be further explored in the future when the prevalence of the disease is higher and when further validation can be done. Strategizing maneuvers, like spaced repeated self-testing, can be done in future studies considering the existing burden of disease.

CONCLUSION

Snort-spit saliva RT-PCR can be an acceptable alternative to nasopharyngeal RT-PCR, especially in a setting where nasopharyngeal specimen cannot be obtained. There is some evidence to suggest that snort-spit saliva RT-PCR may be a more sensitive method than nasopharyngeal RT-PCR. Snort-spit rapid antigen testing is not an acceptable alternative to nasopharyngeal RT-PCR. Further investigation into the superiority of snort-spit saliva RT-PCR over nasopharyngeal RT-PCR is recommended.

Statement of Authorship

All authors participated in the data collection and analysis and approved the final version submitted.

Author Disclosure

All authors declared no conflicts of interest.

Funding Source

This paper was funded by the UP-PGH 2020 COVID Research Grant for Medical Specialist/Consultants.

REFERENCES

1. World Health Organization. Diagnostic testing for SARS-CoV-2. Interim guidance [Internet]. 2020 [cited 2021 Jan]. Available from: <https://www.who.int/publications/i/item/diagnostic-testing-for-sars-cov-2>
2. Xu J, Li Y, Gan F, Du Y, Yao Y. Salivary Glands: Potential Reservoirs for COVID-19 Asymptomatic Infection. *J Dent Res*. 2020 Jul;99(8):989. doi: 10.1177/0022034520918518. Epub 2020 Apr 9. PMID: 32271653.
3. Khurshid Z, Asiri FYI, Al Wadaani H. Human Saliva: Non-Invasive Fluid for Detecting Novel Coronavirus (2019-nCoV). *Int J Environ Res Public Health*. 2020 Mar 26;17(7):2225. doi: 10.3390/ijerph17072225. PMID: 32224986; PMCID: PMC7178089.
4. To KK, Tsang OT, Yip CC, Chan KH, Wu TC, Chan JM, et al. Consistent Detection of 2019 Novel Coronavirus in Saliva. *Clin Infect Dis*. 2020 Jul 28;71(15):841-843. doi: 10.1093/cid/ciaa149. PMID: 32047895; PMCID: PMC7108139.
5. Wylie AL, Fournier J, Casanovas-Massana A, Campbell M, Tokuyama M, Vijayakumar P, et al. Saliva or Nasopharyngeal Swab Specimens for Detection of SARS-CoV-2. *N Engl J Med*. 2020 Sep 24;383(13):1283-1286. doi: 10.1056/NEJMc2016359. Epub 2020 Aug 28. PMID: 32857487; PMCID: PMC7484747.
6. Czumbel LM, Kiss S, Farkas N, Mandel I, Hegyi A, Nagy Á, et al. Saliva as a Candidate for COVID-19 Diagnostic Testing: A Meta-Analysis. *Front Med (Lausanne)*. 2020 Aug 4;7:465. doi: 10.3389/fmed.2020.00465. PMID: 32903849; PMCID: PMC7438940.
7. European Centre for Disease Prevention and Control. Options for the use of rapid antigen tests for COVID-19 in the EU/EEA and the UK [Internet]. 2020 [cited 2021 Jan]. Available from: <https://www.ecdc.europa.eu/en/publications-data/options-use-rapid-antigen-tests-covid-19-eueea-and-uk>
8. RapiGen. BIOCREDIT COVID-19 Ag. RapiGen [Internet]. 2020 [cited 2021 Jan]. Available from: https://www.biovendor.com/file/13419/1-2%20IFU_Covid-19_Ag_Antigen_Multi%20Language%20SEP2020Watermarked.pdf?version=202009300841
9. Bossuyt P, Reitsma J, Bruns D, Gatsonis C, Glasziou P, Irwig L, et al., J. STARD 2015: An updated list of essential items for reporting diagnostic accuracy studies. *BMJ*. 2015 Oct; h5527. <https://doi.org/10.1136/bmj.h5527>
10. Pasomsub E, Watcharananan SP, Boonyawat K, Janchompoo P, Wongtabtim G, Suksuwan W, et al. Saliva sample as a non-invasive specimen for the diagnosis of coronavirus disease 2019: a cross-sectional study. *Clin Microbiol Infect*. 2020 May 15:S1198-743X(20)30278-0. doi: 10.1016/j.cmi.2020.05.001. Epub ahead of print. PMID: 32422408; PMCID: PMC7227531.
11. Altawalrah H, AlHuraish F, Alkandari WA, Ezzikouri S. Saliva specimens for detection of severe acute respiratory syndrome coronavirus 2 in Kuwait: A cross-sectional study. *J Clin Virol*. 2020 Nov;132:104652. doi: 10.1016/j.jcv.2020.104652. Epub 2020 Oct 1. PMID: 33053493; PMCID: PMC7527795.